Trinitron®

Color Video Monitor

Instructions for Use Page 2

Before operating this unit, please read this manual thoroughly and retain it for future reference.

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PVM-20M7MDE

WARNING

To prevent fire or shock hazard, do not expose the unit to rain or moisture.

Dangerously high voltages are present inside the unit.

Do not open the cabinet. Refer servicing to qualified personnel only.

In the event of a malfunction or when maintenance is necessary, consult an authorized Sony dealer.

Power Switch

The power switch is a functional switch only. To isolate the set from the mains supply remove the mains plug from the wall socket.

For customers in the United Kingdom**WARNING**

THIS APPARATUS MUST BE EARTHED

IMPORTANT

The wires in this mains lead are coloured in accordance with the following code:

| | |
|-------------------|---------|
| Green-and-yellow: | Earth |
| Blue: | Neutral |
| Brown: | Live |

As the colours of the wires in the mains lead of this apparatus may not correspond with the coloured markings identifying the terminals in your plug proceed as follows:

The wire which is coloured green-and-yellow must be connected to the terminal in the plug which is marked by the letter E or by the safety earth symbol \perp or coloured green or green-and-yellow.

The wire which is coloured blue must be connected to the terminal which is marked with the letter N or coloured black.


The wire which is coloured brown must be connected to the terminal which is marked with the letter L or coloured red.

Ensure that your equipment is connected correctly - If you are in any doubt consult a qualified electrician.

Important safeguards/notices for use in the medical environments

1. All the equipments connected to this unit shall be certified according to Standard IEC601-1, IEC950, IEC65 or other IEC/ISO Standards applicable to the equipments.
2. When this unit is used together with other equipment in the patient area*, the equipment shall be either powered by an isolation transformer or connected via an additional protective earth terminal to system ground unless it is certified according to Standard IEC601-1 and IEC601-1-1.

* Patient Area

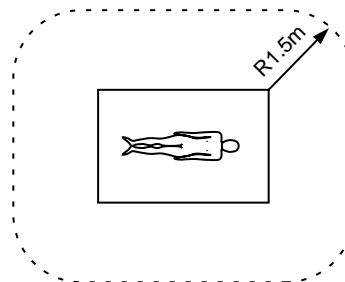3. The leakage current could increase when connected to other equipment.
4. The operator should take precautions to avoid touching the rear panel input and output circuitry and the patient at the same time.
5. Model PVM-20M7MDE is a video monitor intended for use in a medical environment to display video pictures from cameras or other video system.

This unit contains substances which can pollute the environment if disposed carelessly. Please contact our nearest representative office or your local environmental office in case of disposal of this unit.

ATTENTION:

Picture distortion may occur if this monitor is positioned in close proximity to any equipment emitting electromagnetic radiation.

| | |
|--|-----------|
| Precautions | 5 |
| Features | 6 |
| Location and Function of Parts and Controls | 8 |
| Front | 8 |
| Front Operation Panel | 9 |
| Rear Panel | 11 |
| Preparation | 14 |
| Connecting the AC Power Cord | 14 |
| Connecting a Cable to the BNC Connector | 15 |
| Connecting a Cable to the RS-232C Connector | 15 |
| Attaching the Side Covers | 16 |
| Attaching the Control Panel Cover | 16 |
| Operation..... | 17 |
| Getting Started | 17 |
| Adjustment | 18 |
| On-Screen Menus | 19 |
| Using the On-Screen Menus | 19 |
| On-Screen Menu Configuration | 20 |
| Main Menu | 20 |
| USER SERVICE Menus | 24 |
| Troubleshooting | 28 |
| Specifications | 29 |

Symbols on the unit

| Symbol | Location | This symbol indicates |
|---|-----------------|---|
|  | Front panel | Main power switch. Press to turn the monitor on or off. |
|  | Front panel | |
|  | Front panel | |
|  | Rear panel | The equipotential terminal which brings the various parts of a system to the same potential. |
|  | Rear panel | Functional earth terminal |
|  | Rear panel | Alternating current |
|  | Inside the unit | Presence of uninsulated "dangerous voltage" within the product's enclosure that may be sufficient to constitute a risk of electric shock. |
|  | Rear panel | Attention, consult ACCOMPANYING DOCUMENTS |

Attention – when the product is installed in a rack:

- **Elevated operating ambient temperature**
If installed in a closed or multi-unit rack assembly, the operating ambient temperature of the rack environment may be greater than room ambient.
Therefore, consideration should be given to installing the equipment in an environment compatible with the manufacturer's maximum rated ambient temperature of 0 to +40 °C (Tmra).
- **Reduced air flow**
Installation of the equipment in a rack should be such that the amount of air flow required for safe operation of the equipment is not compromised.
- **Mechanical loading**
Mounting of the equipment in the rack should be such that a hazardous condition is not achieved due to uneven mechanical loading.
- **Circuit overloading**
Consideration should be given to the connection of the equipment to the supply circuit and the effect that overloading of circuits might have on overcurrent protection and supply wiring.
Appropriate consideration of equipment nameplate ratings should be used when addressing this concern.
- **Reliable earthing**
Reliable earthing of rack-mounted equipment should be maintained. Particular attention should be given to supply connections other than direct connections to the branch circuit (e.g., use of power strips).

Picture

HR (High Resolution) Trinitron¹⁾ picture tube

The monitor is equipped with the HR Trinitron tube newly developed for professional use in broadcasting stations.

The HR Trinitron tube with Aperture Grill (AG) pitch of 0.30 mm provides a high resolution picture.

Digital comb filter

When PAL video signals are received, a digital comb filter activates to give a more accurate Y/C separation. This contributes to decreasing the vertical cross luminance phenomena which cannot be removed by the analog comb filter. The digital comb filter gives a high resolution picture.

Improvement of field flicker

With the field double scanning system incorporated, the monitor can much reduce flicker phenomena which cause severe eyes fatigue, especially in still image observation and for long time operation.

Input/control connectors

Analog RGB/COMPONENT input connectors

Analog RGB signals or component (Y, R-Y and B-Y) signals from video equipment can be input to these connectors. Input signals can be selected using the operation buttons and the on-screen menu.

External sync input connectors

When an RGB or component signal is input and the sync signal is set to external in the on-screen menu, the monitor can operate according to a sync signal supplied from an external sync generator.

Automatic termination (connector with mark only)

The input connector is terminated at 75 ohms inside when no cable is connected to the loop-through output connector. When a cable is connected to an output connector, the 75-ohm termination is automatically released.

RS-232C interface

The monitor can be remotely controlled by a personal computer via the serial remote interface.

For detailed information on the commands, refer to the supplied interface manual.

Functions

Power saving

On sensing the absence of sync signals for a certain period, the monitor enters stand-by mode to reduce power consumption to approx. 25 %. This contributes to the power saving. You can select either of 30 minutes or 10 minutes from the on-screen menu as a period of sync signals absence.

Automatic reception of a progressive scan (625/50/1:1) signal

The monitor can receive a progressive scan signal when the 625/50/1:1 signal designated in the specifications of the monitor is input. This allows for system expansion.

On-screen menus

Various monitor settings can be changed by using on-screen menus.

Underscan mode

The signal normally scanned outside of the screen can be monitored in the underscan mode. This makes it possible to view an image or screen of data in its entirety.

Auto/manual degauss

Degaussing of the screen is performed automatically when the power is turned on. When the color display becomes non-uniform, degauss the monitor manually by pressing the DEGAUSS button.

Five menu languages

You can select the language used for the on-screen menus from five languages: English, German, French, Italian and Spanish.

1) "Trinitron" is a registered trademark of Sony Corporation.

Side cover and control panel cover

The monitor can be covered with side covers and a control panel cover. The side covers protect the ventilation holes from splashes from medicines and other liquids, and a control panel cover protects the control buttons on the front panel from inadvertent operation.


Quick reference card

The quick reference card is supplied to help you understand the menu configuration and operating method without reading the instruction manual. The quick reference card can be attached to any convenient place with a supplied double-sided adhesive tape.

EIA standard 19-inch rack mounting


Using the SLR-103A slide rail (not supplied), the monitor can be mounted in an EIA standard 19-inch EIA standard rack.

For details on mounting, refer to the instruction manual supplied with the slide rail kit.

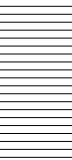Location and Function of Parts and Controls

Front

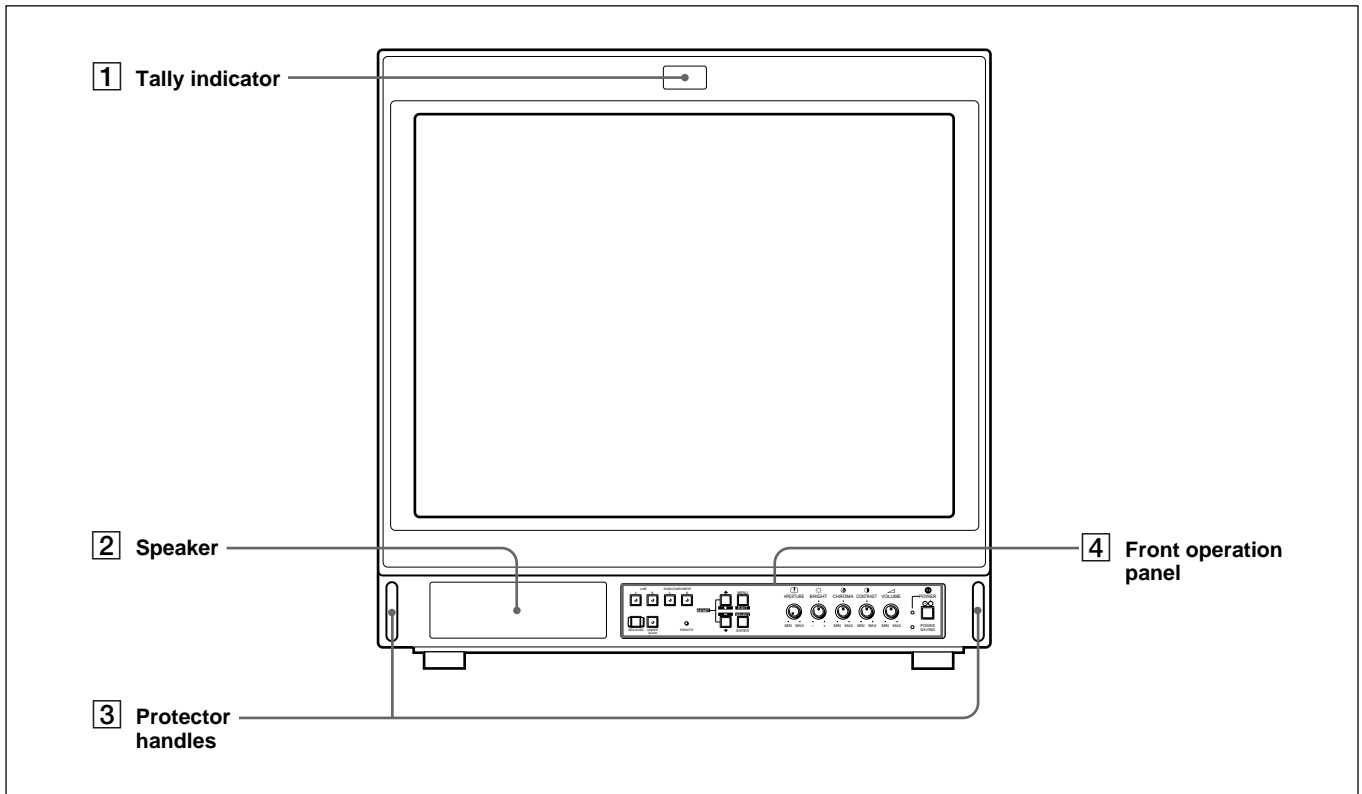1 Tally indicator

Lights when the video equipment connected to this monitor is selected. In order for this indicator to light when the video equipment connected to this monitor is selected, a tally control connection is needed on the RS-232C connector. The light can also be turned on by sending a command from the computer.

For detailed information, see “RS-232C connector (D-sub 9-pin)” in “Specification” on page 31, or refer to the supplied interface manual.

2 Speaker

Provides monaural sounds.

3 Protector handles


Pull out or insert the monitor as holding these handles when mounting the monitor in the EIA standard 19-inch rack. These handles also function to protect the control buttons and knobs on the front panel from being broken when you tilts the monitor.

4 Front operation panel

Operation buttons, control knobs and indicators are located on the front operation panel.

For detailed information on the operation panel, see “Front Operation Panel” on the next page.

Front Operation Panel

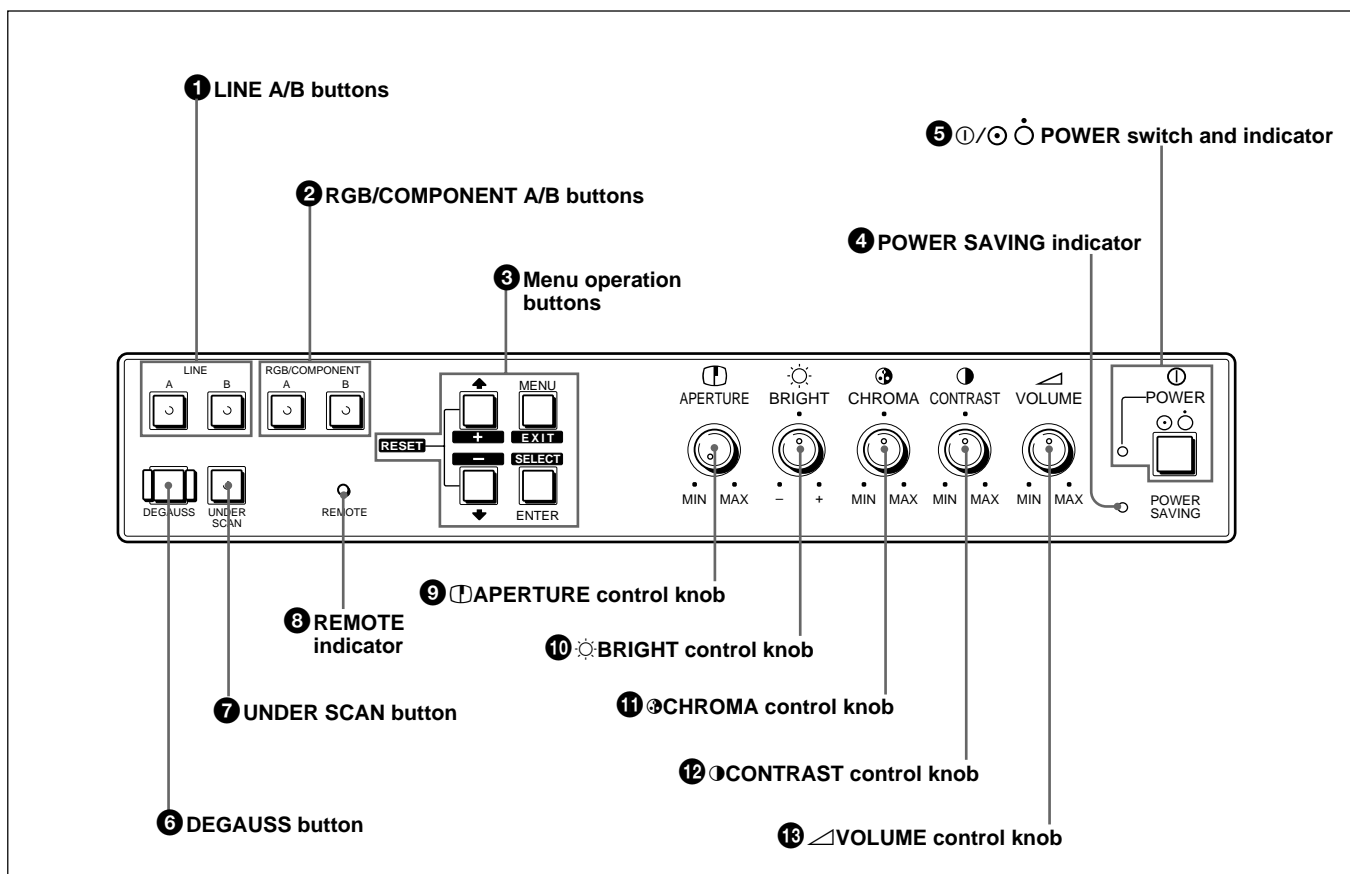① LINE A/B (input select) buttons

Press to select a signal to be monitored. The indicator on the pressed button lights.

| Input signal | Press |
|---|--------|
| Signal fed through the LINE A connector | LINE A |
| Signal fed through the LINE B connector | LINE B |

② RGB/COMPONENT A/B (input select) buttons

Press to select a signal to be monitored. The indicator on the pressed button lights.

| Input signal | Press |
|---|-----------------|
| Signal fed through the RGB/COMPONENT A connectors ^{a)} | RGB/COMPONENT A |
| Signal fed through the RGB/COMPONENT B connectors ^{a)} | RGB/COMPONENT B |

a) The signal fed through the RGB/COMPONENT A/B connectors can be selected from the RGB SYSTEM menu.

For detailed information on selecting the RGB or COMPONENT signal, see “To select the signal to be monitored through the RGB/COMPONENT connectors” on page22.


Location and Function of Parts and Controls

③ Menu operation buttons

Press to make the menu appear.

To return to the previous screen, press the MENU button again.

For detailed information on menu operation buttons, see “Menu operation buttons” on page 19.

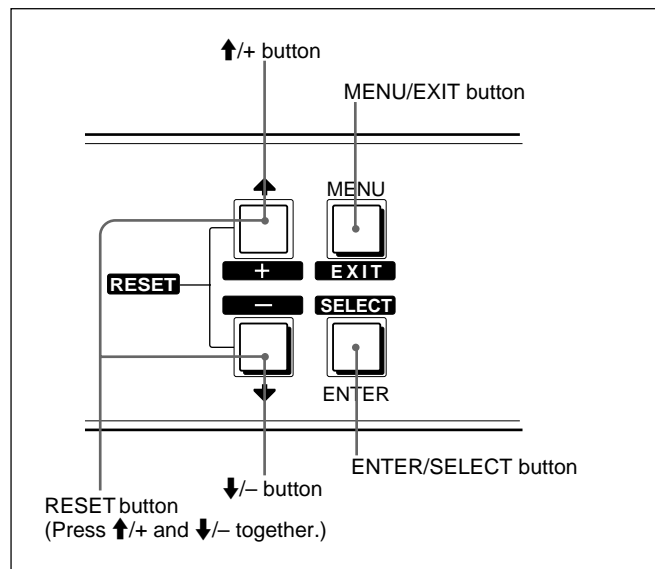④ POWER SAVING indicator

This indicator lights when the POWER SAVING function is active.

⑤ ①/② Ⓞ POWER switch and indicator

Press to turn the monitor on. The green power indicator lights.

Press again to turn the power off.

⑥ DEGAUSS button

Press this button momentarily to demagnetize the screen.

Wait for at least 10 minutes before pressing this button again.

⑦ UNDER SCAN button

Press this button to display the whole picture. The button indicator lights. The display size is reduced by approximately 5 % so that the four corners of the raster are visible.

Press this button again to reset the display to the normal size. The button indicator goes out.

⑧ REMOTE indicator

This indicator lights when REMOTE (RS-232C) is set to REMOTE or REMOTE & LOCAL on the menu.

⑨ ① APERTURE control knob

Adjusts the sharpness of the screen.

For how to make adjustment, see “Adjustment” on page 18.

⑩ ① BRIGHT (brightness) control knob

Adjusts the brightness of the screen.

For how to make adjustment, see “Adjustment” on page 18.

⑪ ① CHROMA (chrominance) control knob

Adjusts the color intensity of the video signal.

For how to make adjustment, see “Adjustment” on page 18.

⑫ ① CONTRAST control knob

Adjusts the contrast of the screen.


For how to make adjustment, see “Adjustment” on page 18.

⑬ ① VOLUME control knob

Adjusts the speaker volume.

For how to make adjustment, see “Adjustment” on page 18.

Rear Panel

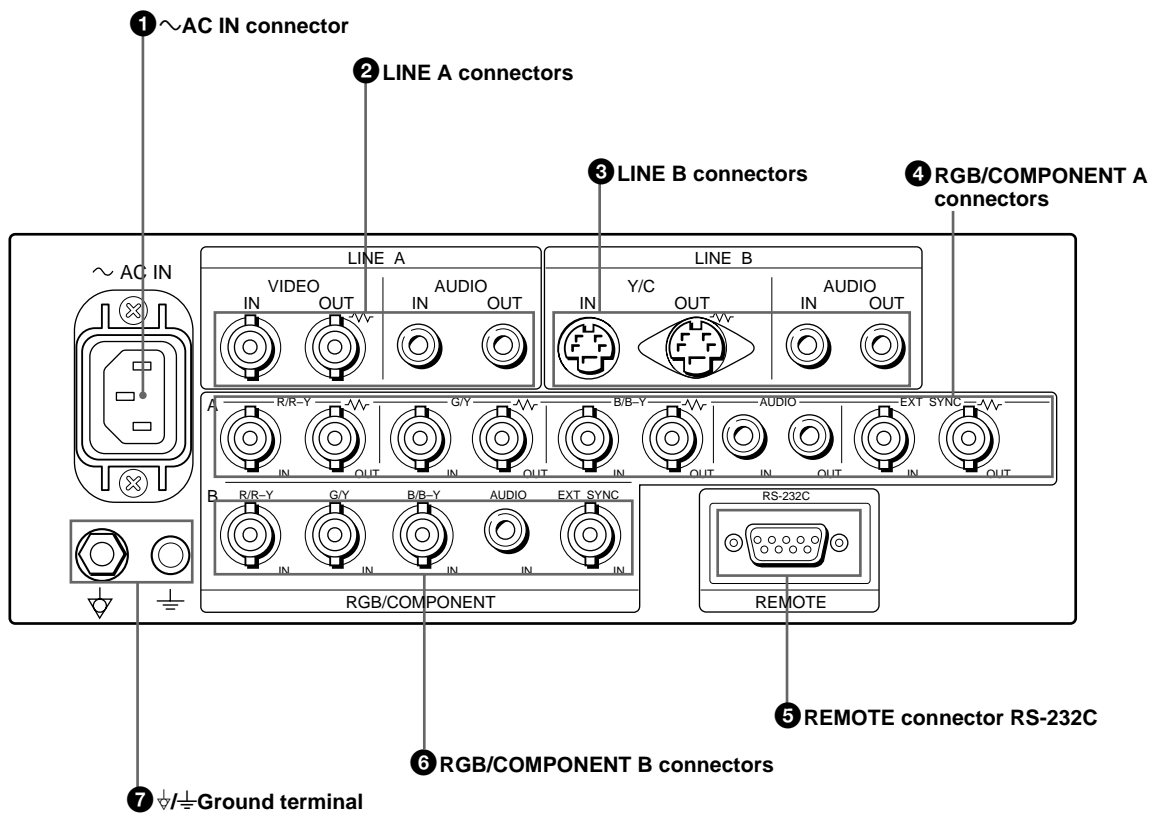Note

Before connecting the video equipment, see "Important safeguards/ notices for use in the medical environments" on page 2.

(The \sim mark indicates automatic termination.)

1 ~AC IN (inlet) connector

Connect the supplied AC power cord to this connector and to a wall outlet.

2 LINE A connectors

Input connectors for the composite video and audio signals and their loop-through output connectors. To monitor the input signal fed to these connectors, press the LINE A button on the front panel.

VIDEO IN (input) connector (BNC-type)

Connect to the video output connector on video equipment, such as a VCR or a color video camera. For a loop-through connection, connect to the video output connector on another monitor.

VIDEO OUT (output) connector (BNC-type)

Loop-through output connector of the VIDEO IN connector. Connect to the video input connector on a VCR or another monitor.

When the cable is connected to this connector, the 75-ohm input termination is automatically released, and the signal input to the VIDEO IN connector is output from this connector.

AUDIO IN (input) connector (phono jack)

Connect to the audio output connector on a VCR or other equipment. For a loop-through connection, connect to the audio output on another monitor.

AUDIO OUT (output) connector (phono jack)

Loop-through output on the AUDIO IN connector. Connect to the audio input connector of a VCR or another monitor.

Location and Function of Parts and Controls

③ LINE B connectors

Input connectors for the Y/C separate video and audio signals and their loop-through output connectors.

To monitor the input signal fed through these connectors, press the LINE B button on the front panel.

Y/C IN (input) connector (4-pin mini-DIN)

Connect to the Y/C separate output connector of a color video camera, VCR or other video equipment.

For pin assignments of the Y/C IN connector, see “Y/C IN connector (4-pin mini-DIN)” on page 30.

Y/C OUT (output) connector (4-pin mini-DIN)

Loop-through output of the Y/C IN connector.

Connect to the Y/C separate input connector on a VCR or another monitor.

When the cable is connected to this connector, the 75-ohm input termination is automatically released, and the signal input to the Y/C IN connector is output from this connector.

AUDIO IN (input) connector (phono jack)

Connect to the audio output connector on a VCR or other equipment. For a loop-through connection, connect to the audio output on another monitor.

AUDIO OUT (output) connector (phono jack)

Loop-through output of the AUDIO IN connector. Connect to the audio input connector on a VCR or another monitor.

④ RGB/COMPONENT A connectors

Input connectors for the R/G/B signals or component signals, external sync signals and audio signals and their loop-through output connectors.

To monitor the input signal fed through these connectors, press the RGB/COMPONENT A button on the front panel. Then, select one of the four items RGB -EXT SYNC, RGB -SYNC ON G, COMP -EXT SYNC or COMP -SYNC ON Y on the RGB SYSTEM -RGB A- menu.

For detailed information on settings, see “To select the signal to be monitored through the RGB/COMPONENT connectors” on page 22.

R/G/B IN or R-Y/Y/B-Y IN (input) connectors (BNC-type)

When RGB -SYNC ON G or COMP -SYNC ON Y is selected on the RGB SYSTEM -RGB A- menu, the monitor operates on the sync signal from the G/Y channel.

When RGB-EXT SYNC or COMP-EXT SYNC is selected on the RGB SYSTEM -RGB A- menu, the monitor operates on an external sync signal.

To monitor the analog RGB signal: Connect to the analog RGB output connectors on a video camera, VCR or other video equipment.

To monitor the component signal: Connect to the R-Y/Y/B-Y component output connectors on video equipment.

R/G/B OUT or R-Y/Y/B-Y OUT (output) connectors (BNC-type)

Loop-through outputs of the R/R-Y, G/Y IN and B/B-Y IN connectors.

When cables are connected to these connectors, the 75-ohm input termination is automatically released, and the signal inputs to the R/R-Y, G/Y IN and B/B-Y IN connectors are output from these connectors.

To output the analog RGB signal: Connect to the analog RGB input connectors on a video printer or another monitor.

To output the component signal: Connect to the R-Y/Y/B-Y component input connectors on video equipment.

AUDIO IN (input) connector (phono jack)

Connect to the audio output connectors on video equipment when the analog RGB or component signal is input.

AUDIO OUT (output) connector (phono jack)

Loop-through output of the AUDIO IN connector. Connect to the audio input connector on a VCR or another monitor.

EXT SYNC IN (external sync input) connector (BNC-type)

Connect to the sync signal output on a video camera, VCR or other video equipment, when this monitor operates on an external sync signal.

To use the sync signal fed through this connector, select RGB -EXT SYNC or COMP -EXT SYNC on the RGB SYSTEM -RGB A- menu.

EXT SYNC OUT (external sync output) connector (BNC-type)

Loop-through output of the EXT SYNC IN connector. Connect to the external sync input connector on video equipment to be synchronized with this monitor.

When the cable is connected to this connector, the 75-ohm input termination is released, and the signal input to the EXT SYNC IN connector is output from this connector.

For detailed information on settings, see “To select the signal to be monitored through the RGB/COMPONENT connectors” on page 22.

⑤ REMOTE connector RS-232C (D-sub 9-pin)

Connect to the RS-232C control connector on a personal computer. You can operate the monitor with the commands from the computer.

For detailed information, refer to the supplied “Interface Manual for Programmers”.

For pin assignments of the REMOTE connector, see “RS-232C connector (D-sub 9-pin)” on page 31.

⑥ RGB/COMPONENT B connectors

Input connectors for the R/G/B signals or component signals, external sync signals and audio signals.

To monitor the input signal fed through these connectors, press the RGB/COMPONENT B button on the front panel. Then, select one of the four items RGB -EXT SYNC, RGB -SYNC ON G, COMP -EXT SYNC or COMP -SYNC ON Y on the RGB SYSTEM -RGB B- menu.

For detailed information on settings, see “To select the signal to be monitored through the RGB/COMPONENT connectors” on page 22.

R/G/B IN or R-Y/Y/B-Y IN (input) connectors (BNC-type)

When RGB -SYNC ON G or COMP -SYNC ON Y is selected on the RGB SYSTEM -RGB B- menu, the monitor operates on a sync signal from the G/Y channel.

When RGB -EXT SYNC or COMP -EXT SYNC is selected on the RGB SYSTEM -RGB B- menu, the monitor operates on the external sync signal.

To monitor the analog RGB signal: Connect to the analog RGB output connectors on a video camera, VCR or other video equipment.

To monitor the component signal: Connect to the R-Y/Y/B-Y component output connectors on video equipment.

AUDIO IN (input) connector (phono jack)

Connect to the audio output connectors on video equipment when the analog RGB or component signal is input.

EXT SYNC IN (external sync input) connector (BNC-type)

Connect to the sync signal output on a video camera, VCR or other video equipment, when this monitor operates on an external sync signal.

To use the sync signal fed through this connector, select RGB -EXT SYNC or COMP -EXT SYNC on the RGB SYSTEM -RGB B- menu.

⑦ Ground terminal

Equipotential terminal

Connect to the equipotential plug to bring the various parts of a system to the same potential.

Functional earth

See “Symbols on the unit” on page 4.

Preparation

This section explains preparation required before attempting to operate the monitor, and how to connect to other equipment.

For detailed information on connecting the monitor to other equipment, refer also to the instruction manuals supplied with those equipment.

Note

Always turn off the power of the monitor and other equipment before attempting to make connections.

Connecting the AC Power Cord


Connect the supplied AC power cord to the \sim AC IN connector and to a wall outlet.

Cautions on power connection

- Use the supplied AC power cord.
- Set the power switch to OFF before connecting or removing the AC power cord.
- Make sure that the power supply conforms to the voltage rating in “Specifications” on page 29.

To connect an AC power cord securely with the AC plug holder

Before connecting the AC power cord to the \sim AC IN connector and to a wall outlet, we recommend to connect an AC power cord securely with the AC plug holder as follows .

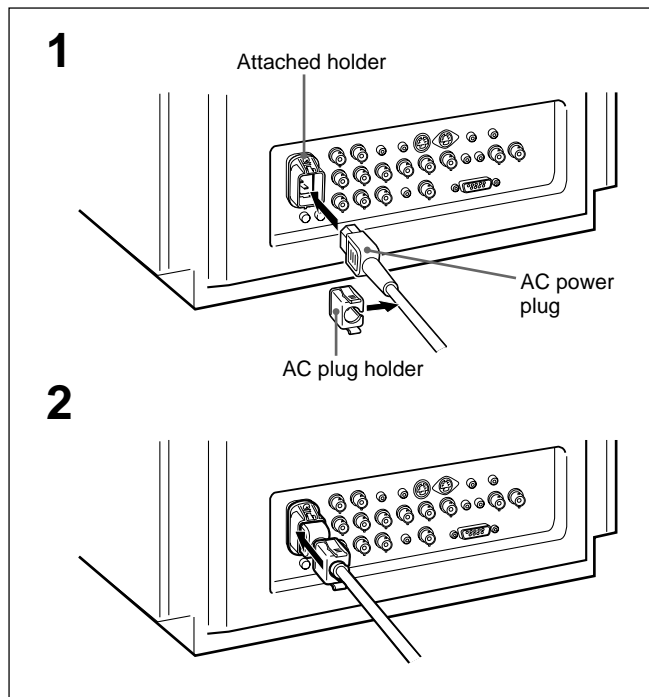- 1** Plug the AC power plug into the attached holder. Then attach the supplied AC plug holder on the top of the AC power cord.
- 2** Slide the AC plug holder over the cord until it connects with the attached holder.


To disconnect the AC power cord

Note

Wait at least 10 seconds after switching OFF the power switch before disconnecting the AC power cord to discharge any static electricity from the CRT display tube.

Pull out AC plug holder by squeezing the up and down side.

Connecting a Cable to the BNC Connector

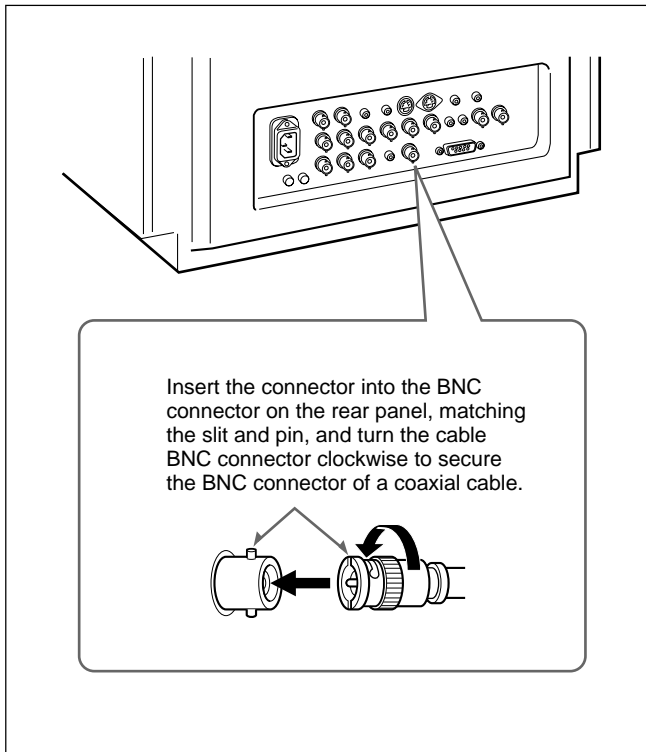
Note

Before connecting the video equipment, see **“Important safeguards/notices for use in the medical environments”** on page 2.

Connect a coaxial cable with the BNC connectors to the BNC connectors on the rear panel as illustrated.

Connecting a Cable to the RS-232C Connector

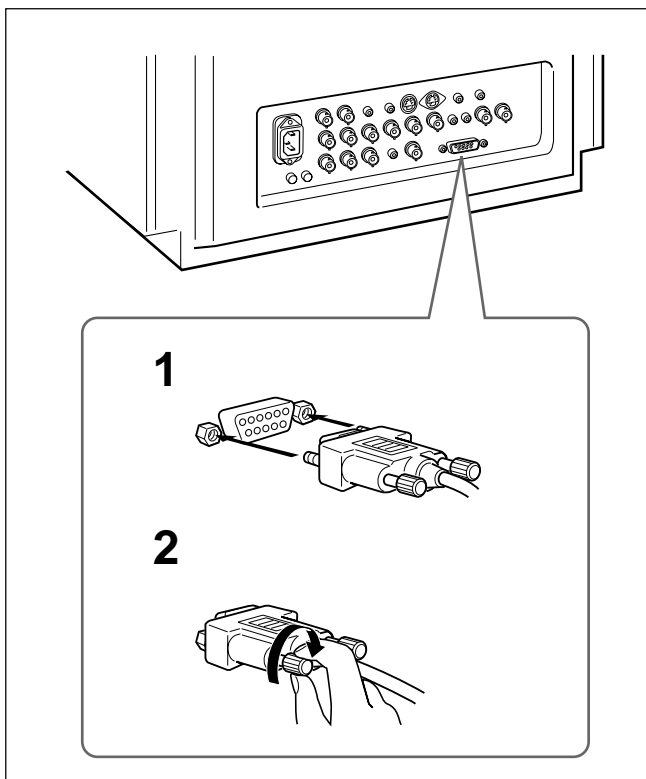Note


Before connecting the video equipment, see **“Important safeguards/notices for use in the medical environments”** on page 2.

- 1 Align the plug and socket carefully.
- 2 Push fully in and tighten the screws by hand.

To disconnect the plug

Loosen the screws, and pull out the plug.

Attaching the Side Covers

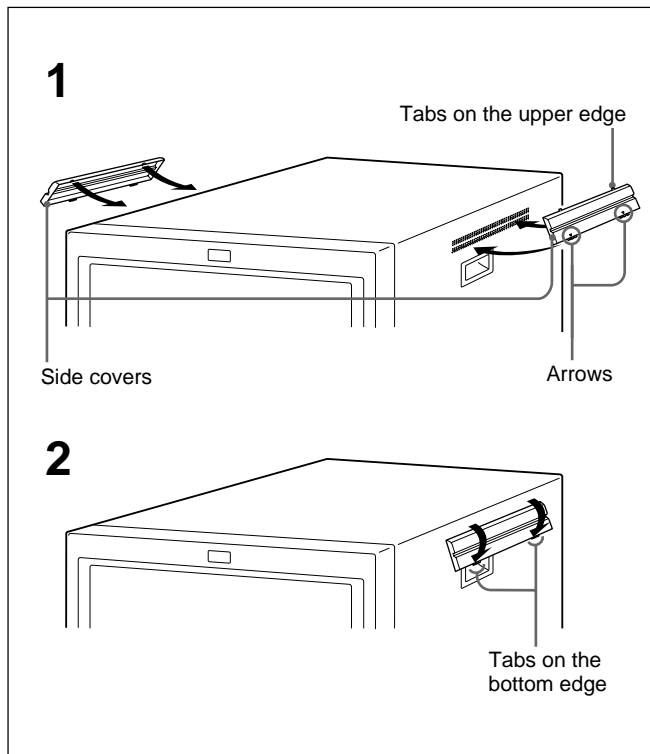In order to protect the ventilation holes from splashes from medicines, etc., attach the supplied side covers as illustrated.

- 1 Hook the tabs on the upper edge into the ventilation holes, making sure that the arrows on the cover are facing down.


Note

Attach the side covers on all ventilation holes.

- 2 Push up the tabs on the bottom edge and fit the cover into the lowest ventilation holes.

Attach covers on both left and right vents.


Attaching the Control Panel Cover

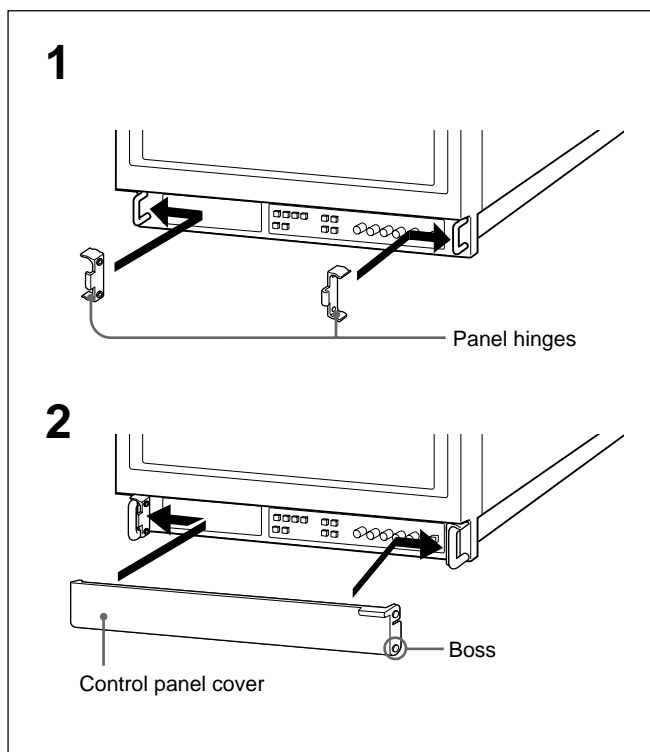In order to protect the control buttons on the front panel from inadvertent operation, attach the supplied control panel cover.

- 1 Attach the panel hinges to the left and right grips from the inner side.
- 2 Fit the bosses on both sides of the control panel cover into the lower holes on the panel hinges, bending the control panel cover a little.

Getting Started

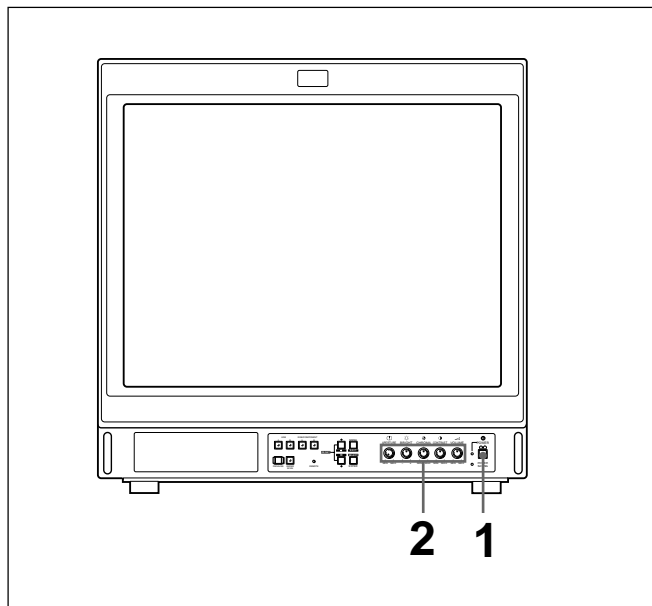Notes

- After the power has been turned on or you have pressed the DEGAUSS button on the front panel of the monitor, the monitor is demagnetized for approximately 10 seconds. This generates a strong magnetic field around the monitor which may affect data stored on magnetic tape or disks near the monitor. Place all magnetic recording equipment and tapes/disks well away from the monitor.
- Whenever you change the location of the monitor, degauss the monitor before adjustments are made.

- 1** Set the power switch to ON and then turn on the power on the other equipment.

We recommend warming up the monitor for at least 30 minutes after turning the power switch on and inputting the signal.

- 2** Adjust the monitor settings as required.

For detailed information on adjustments, see "Adjustment" on page 18, "Main menu functions" on page 22 and "USER SERVICE menu functions" on page 26.

When the POWER SAVING indicator is lit

If no sync signals are received for a certain period, the POWER SAVING indicator lights, and the monitor screen blanks out. This is power saving mode. Touching any button other than the POWER switch will restore the monitor to normal operation.

Note

It takes about ten seconds for the picture to reappear on the screen after power saving mode is released.

If the color display becomes non-uniform

Press the DEGAUSS button so that the color uniformity becomes as even. If the color display becomes non-uniform even after pressing the DEGAUSS button, landing adjustment is required on the on-screen menu.

For details, see "To adjust the landing" on page 27.

To display the whole picture on the monitor

Press the UNDER SCAN button so that the signal normally scanned outside of the screen appears on the monitor.


Note

When the monitor is in underscan mode, dark RGB scanning lines may appear on the top edge of the screen. These are caused by an internal test signal, and are not part of the input signal.

Operation

Adjustment

This section explains adjustment performed using controls on the front panel of the monitor.

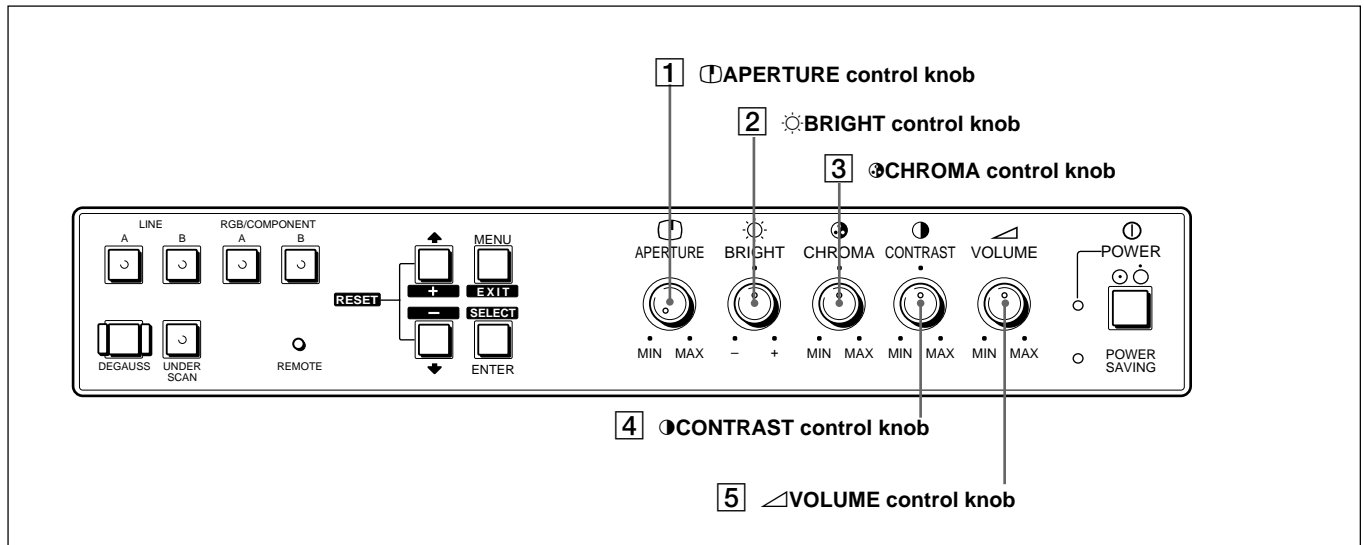To adjust the sharpness

Use the **APERTURE** control **1**.
Turn the control towards MAX to sharpen the picture.
Turn the control towards MIN to soften the picture.

To adjust the brightness

Use the **BRIGHT** control **2**.
Turn the control towards + to make the picture brighter.
Turn the control towards - to make the picture darker.

To adjust the color intensity of the video signal (color saturation)

Use the **CHROMA** control **3**.
Turn the control towards MAX to increase the color intensity.
Turn the control towards MIN to decrease the color intensity.

To adjust the contrast

Use the **CONTRAST** control **4**.
Turn the control towards MAX to increase the contrast.
Turn the control towards MIN to decrease the contrast.

To adjust the speaker volume

Use the **VOLUME** control **5**.
Turn the control towards MAX to increase the volume.
Turn the control towards MIN to decrease the volume.


On-Screen Menus

The on-screen menus allow you to make various settings and adjustments of the monitor. The settings remain in memory even when the power is turned off.

Using the On-Screen Menus

Menu operation buttons

There are four menu operation buttons on the front panel of the monitor.

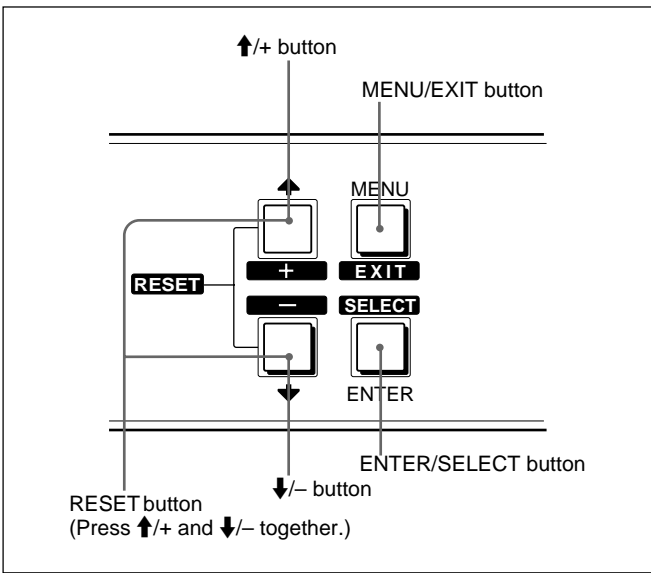To display the on-screen menus

To display the main menu (called MENU on the on-screen menu), press the MENU button. To display the USER SERVICE menu, press and hold down the MENU button for 2 or 3 seconds.

To exit the on-screen menus

Each time you press the MENU/EXIT button, the screen returns to the one previously displayed. Press the MENU button until the regular screen appears.


Functions of menu buttons

The function of the menu operation buttons depends on the purpose that you want to do on the displayed on-screen menu.

The following table shows the button functions.

| Button | To select menu item | To adjust selected menu item |
|--------------------------------------|--|--|
| MENU EXIT | Exits the menu and returns to the previous menu. | Exits the menu and returns to the previous menu. |
| ENTER SELECT | Decides a selected item. | Select an item, only on the ADJUST GAIN and ADJUST BIAS MENU menu. |
| ↑ + | Moves the cursor (▶) upwards. | Increases the selected value. |
| ↓ - | Moves the cursor (▶) downwards. | Decreases the selected value. |
| RESET Press ↑/+ and ↓/- together. | | Can be used only on ADJUST GAIN, ADJUST BIAS and LANDING menus. Resets the setting to the original value before this change on ADJUST GAIN and ADJUST BIAS menus. Resets the setting to the standard value (00) on the LANDING menu. |

Buttons that can be used on each menu are displayed on the bottom line of the screen.

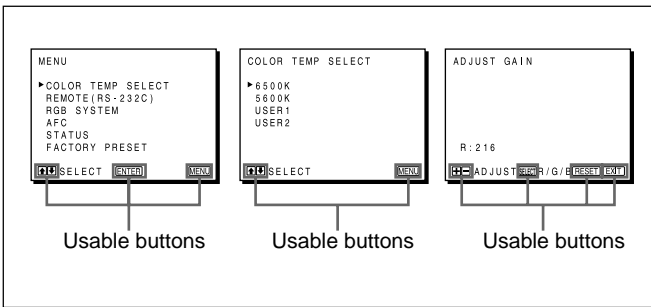Usable menu button display line

On-Screen Menus

On-Screen Menu Configuration


The monitor has two menus: Main menu (called MENU on the on-screen menu) and USER SERVICE menu.

These two menus have sub menus.

For detailed information on the configuration and functions of these on-screen menus, see the following “Main Menu” and “USER SERVICE Menus” on page 24.

Main Menu

Main menu tree

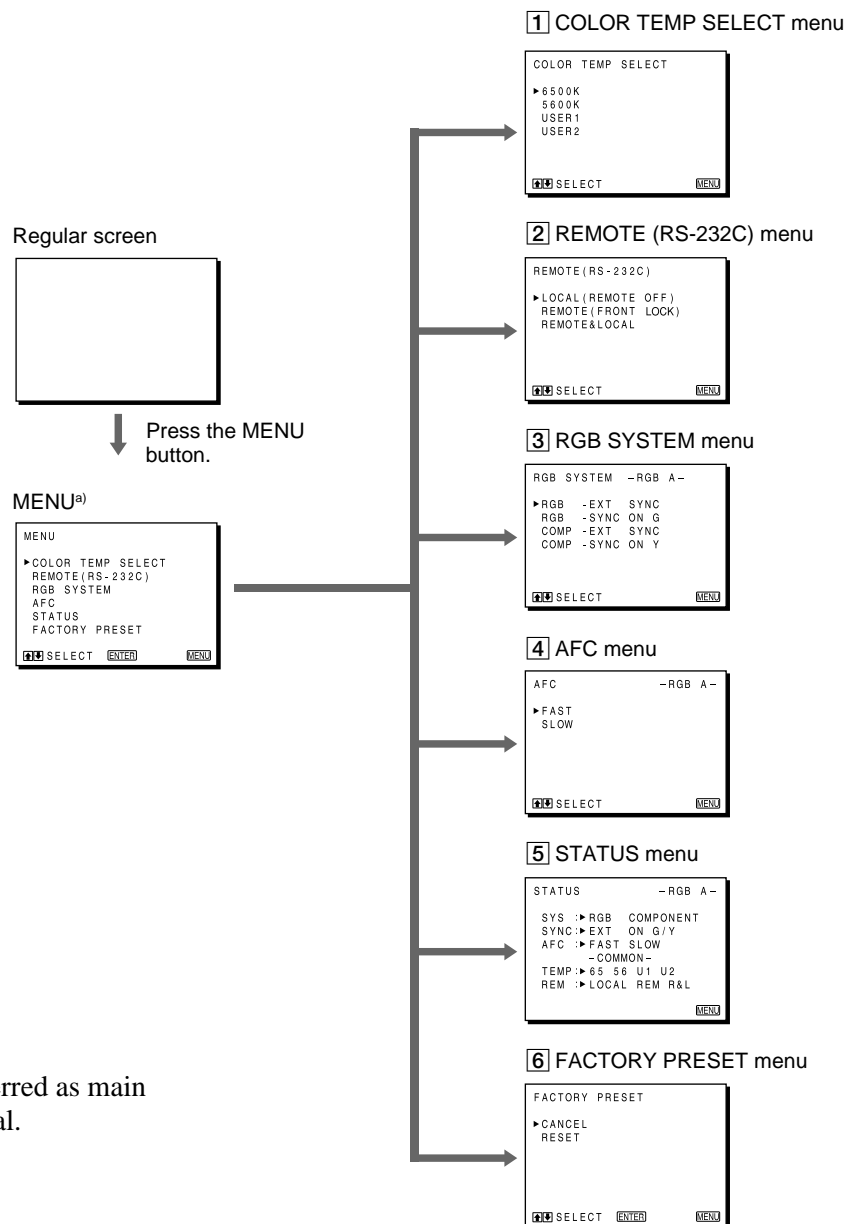List of operation menus

The following list shows the function and factory setting and reference page of operation menus. You can change the factory setting using the menu. To reset the settings to the factory-settings, select RESET from the FACTORY PRESET menu.

For details, see “To reset the monitor status to the factory settings” on page 23.

| Operation menu | Function | Factory setting | Reference page |
|-----------------------|--|--------------------|----------------|
| [1] COLOR TEMP SELECT | Selects the color temperature | 6500K | 22 |
| [2] REMOTE (RS-232C) | Selects monitor operation mode when controlling the monitor via the RS-232C interface | LOCAL (REMOTE OFF) | 22 |
| [3] RGB SYSTEM | Selects the signal to be monitored through the RGB/COMPONENT connectors | RGB-EXT SYNC | 22 |
| [4] AFC | Selects the AFC time constant to correct the picture of VCR when it is skew on the upper part of the monitor screen. | FAST | 23 |
| [5] STATUS | Checks the current menu settings for the monitor | | 23 |
| [6] FACTORY PRESET | Resets the monitor status to the factory settings. | CANCEL | 23 |

On-Screen Menus

Displaying the main menu

Press the MENU button.
The main menu appears.


Main menu functions

This section explains settings which you can change using the operation menus.

To select the color temperature

The monitor has two preset color temperatures, 6500K and 5600K and two customized color temperatures USER 1 and USER 2. You can select one of those four standard temperatures on the COLOR TEMP SELECT menu.

1 COLOR TEMP SELECT menu
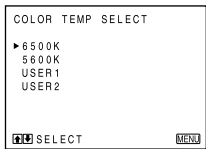
| Item | Function |
|--------|--|
| 6500K | Sets the color temperature to 6500K. |
| 5600K | Sets the color temperature to 5600K. |
| USER 1 | Sets the color temperature to the one set on USER 1 COLOR TEMP ADJ menu. |
| USER 2 | Sets the color temperature to the one set on USER 2 COLOR TEMP ADJ menu. |


You can adjust the color temperature on the USER SERVICE menus of USER 1 and USER 2 to obtain a color temperature other than the factory settings of 6500K and 5600K.

For details of how to obtain a color temperature other than the factory settings, see "To adjust the white balance" on page 26.

To control the monitor via the RS-232C interface

You can select the monitor operation mode on the REMOTE (RS-232C) menu.

2 REMOTE (RS-232C) menu

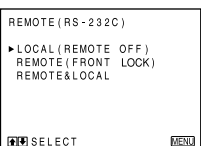
| Item | Function |
|---------------------|--|
| LOCAL (REMOTE OFF) | Operates the monitor using the buttons and control knobs on the front panel of the monitor. |
| REMOTE (FRONT LOCK) | Controls the monitor remotely from a computer. Buttons and control knobs except the POWER switch do not function. To return to the REMOTE (RS-232C) menu from this mode, press the MENU/EXIT button for about 2 seconds. |
| REMOTE & LOCAL | Controls the monitor remotely from a computer. Control knobs such as APERTURE, BRIGHT, CHROMA, CONTRAST and VOLUME do not function. |

To select the signal to be monitored through the RGB/COMPONENT connectors

You can select the signal to be monitored through the RGB/COMPONENT connectors on the RGB SYSTEM menu.

3 RGB SYSTEM menu

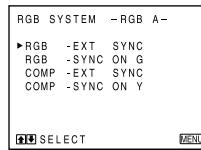
This menu appears when you press the RGB/COMPONENT A button on the front panel. When you press the RGB/COMPONENT B button, the RGB SYSTEM -RGB B- menu appears.

| Item | Function |
|-----------------|--|
| RGB -EXT SYNC | Monitors the RGB signal while operating the monitor according to an external sync signal fed through the EXT SYNC connector. |
| RGB -SYNC ON G | Monitors the RGB signal while operating the monitor according to the sync signal from G channel. |
| COMP -EXT SYNC | Monitors the component signal while operating the monitor according to an external sync signal fed through the EXT SYNC connector. |
| COMP -SYNC ON Y | Monitors the component signal while operating the monitor according to the sync signal from Y channel. |

When the LINE A/B button is pressed


If you select the RGB SYSTEM menu when LINE A or LINE B is selected with the LINE A/B button, the following menu appears.

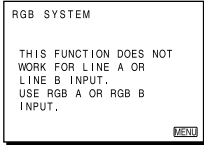To select the AFC (Automatic Frequency Control) time constant

If the signal input to the monitor from a VCR is deviated from the color system specifications, vertical lines of the picture may skew on the upper part of the monitor screen. You can reduce this problem by selecting the proper AFC time constant on the AFC menu.

4 AFC menu
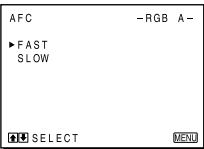
| Item | Function |
|------|---------------------------------------|
| FAST | Sets the AFC time constant to 1 msec. |
| SLOW | Sets the AFC time constant to 2 msec. |


Note

If the picture is still skew even if you change the AFC setting, another problem can be considered. Check your VCR.

To confirm the current menu settings for the monitor

You can confirm the status of the current settings on each menu on the STATUS menu.

5 STATUS menu

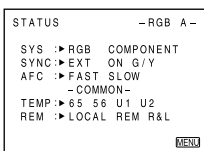When LINE A or LINE B is selected with the LINE A/B button

When the LINE A or LINE B is selected, the column of SYS and SYNC are displayed as follows.


SYS : -- --

SYNC: -- --

To reset the monitor status to the factory settings

You can reset the settings on all menus to the factory settings on the FACTORY PRESET menu.


6 FACTORY PRESET menu
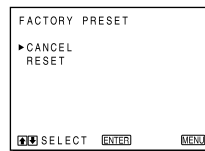
| Item | Function |
|--------|--|
| CANCEL | Cancels "FACTORY PRESET". |
| RESET | Resets all settings to the factory settings except for the "LANGUAGE" setting. |

USER SERVICE Menus

USER SERVICE menu tree

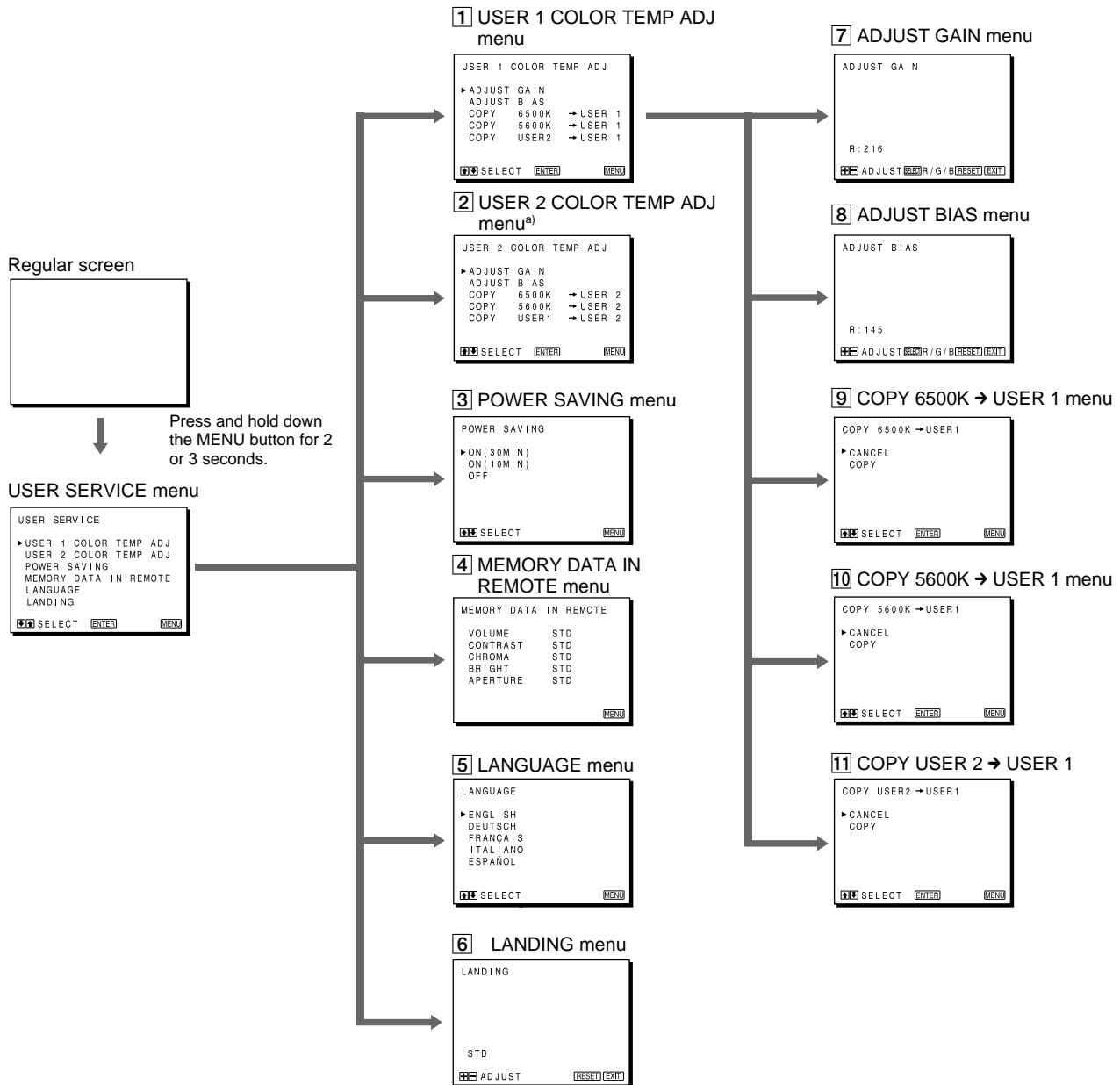a) This menu has the same sub menus as the USER 1 COLOR TEMP ADJUST menu: ADJUST GAIN, ADJUST BIAS, COPY 6500K → USER 2, COPY 5600 K → USER 2 and COPY USER 1 → USER 2.

List of operation menus

The following list shows the function and factory setting and reference page of USER SERVICE menus. You can change the factory setting using the menu. To reset the settings to the factory-settings, select RESET from the FACTORY PRESET menu.

For details, see “To reset the monitor status to the factory settings” on page 23.

| User service menu | Function | Factory setting | Reference page |
|---------------------------|---|------------------------|----------------|
| [1] USER 1 COLOR TEMP ADJ | Adjusts the color temperature of USER 1 setting. | ADJUST GAIN | 26 |
| [2] USER 2 COLOR TEMP ADJ | Adjusts the color temperature of USER 2 setting. | ADJUST GAIN | 26 |
| [3] POWER SAVING | Activates the power saving function. | ON (30 MIN) | 27 |
| [4] MEMORY DATA IN REMOTE | Confirms the data set by the computer in remote control mode. | | 27 |
| [5] LANGUAGE | Selects the menu language. | ENGLISH | 27 |
| [6] LANDING | Adjusts the signal landing. | STD | 27 |
| [7] ADJUST GAIN | Adjusts the gain of R, G or B channel. | USER 1: data of 6500 K | 26 |
| [8] ADJUST BIAS | Adjusts the bias of R, G or B channel. | USER 2: data of 5600 K | 26 |
| [9] COPY 6500K → USER 1 | Copies the color temperature data of 6500 K to USER 1. | CANCEL | 26 |
| [10] COPY 5600K → USER 1 | Copies the color temperature data of 5600 K to USER 1. | CANCEL | 26 |
| [11] COPY USER 2 → USER 1 | Copies the color temperature data of USER 1 to USER 2. | CANCEL | 26 |

On-Screen Menus

Displaying the USER SERVICE menu

With the regular screen displayed, press and hold down the MENU button for 2 or 3 seconds.

The USER SERVICE menu appears.

USER SERVICE menu functions

This section explains settings which you can change using the USER SERVICE menus.

To adjust the white balance

The white balance can be adjusted using the color temperature setting.


You can adjust the color temperature on the USER SERVICE menus of USER 1 and USER 2 to obtain a color temperature other than the factory settings of 6500K and 5600K.

Note

USER 1 is set to 6500K and USER 2 is set to 5600K in the factory setting each.

The following explains how to change the USER 1 COLOR TEMP ADJ menu setting.

1 USER 1 COLOR TEMP ADJ menu

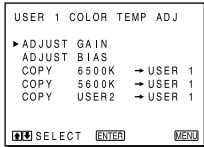
- 1 Select the desired color temperatures (6500K, 5600K or USER 2) to be used as standard when adjusting the color temperature by pressing the **↑** or **↓**.

| Item | Meaning |
|----------------------|--|
| COPY 6500K → USER 1 | The color temperature data of 6500K is copied to USER 1. |
| COPY 5600K → USER 1 | The color temperature data of 5600K is copied to USER 1. |
| COPY USER 2 → USER 1 | The color temperature data of USER 1 is copied to the one set as USER 2. ^{a)} |

a) This is effective to preserve the original set temperature of USER 2. You can modify the value based on the values of USER 2.

- 2 Press the ENTER button.
The following menu appears.

9 COPY 6500K → USER 1 menu

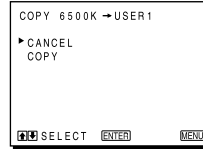- 3 Move the cursor (▶) to COPY by pressing the **↑** or **↓** and press the ENTER button.

The message “DATA COPIED” appears when copying is completed.


- 4 Press the MENU button.
The standard color temperature menu displayed in step 1 appears.

- 5 Adjust the bias and gain if you want to modify the copied color temperature.

| Item | Function |
|------|---|
| GAIN | To adjust white balance (used to compensate for high intensity lighting) |
| BIAS | To adjust the black balance (used to compensate for low intensity lighting) |

- ① Move the cursor (▶) to ADJUST GAIN by pressing the **↑** or **↓** and press the ENTER button.
The ADJUST GAIN menu appears.

7 ADJUST GAIN menu

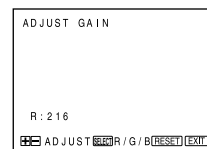- ② Select the R, G or B channel by pressing the SELECT button.
- ③ Adjust the volume of the selected channel by pressing the + or – buttons.
- ④ Press the EXIT button.
The standard color temperature menu reappears and the cursor is positioned on ADJUST BIAS.

- ⑤ Repeat step ② and ③ to adjust the bias.
The screen displays with the adjusted value, but the data is not yet stored.

To reset the settings to the previous values


Press the ↑/+ and ↓/- buttons (these buttons work as RESET button on the menu) simultaneously.

- 6 Press the MENU button.
The adjusted values are stored in memory.

To activate the power saving function

Turn the power saving function on in the following menu.

3 POWER SAVING menu

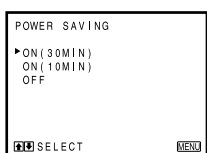
| Item | Function |
|-------------|--|
| ON (30 MIN) | Reduces the power consumption to approx. 25%, 30 minutes after sensing the absence of sync signal. |
| ON (10 MIN) | Reduces the power consumption to approx. 25%, 10 minutes after sensing the absence of sync signal. |
| OFF | Turns off the power saving function. |

To confirm the data set by the computer in remote control mode

You can confirm the settings controlled from the remote computer, such as volume, contrast, chroma, brightness and aperture. When you select MEMORY DATA IN REMOTE on the USER SERVICE menu and press the ENTER button, the MEMORY DATA IN REMOTE menu appears.

4 MEMORY DATA IN REMOTE menu

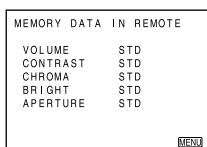To select the menu languages

You can select the menu language from the following five languages: English, German, French, Italian and Spanish on the LANGUAGE menu.


5 LANGUAGE menu

To adjust the landing

If the color is not uniform even after you press the DEGAUSS button, you can adjust the landing so as to obtain color uniformity on the LANDING menu.

6 LANDING menu

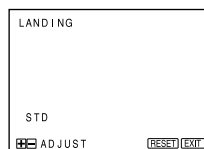The following two methods are available to adjust the landing.

One is to adjust landing as inputting the signal with which the horizontal lines are displayed.

The other is to adjust landing as inputting the signal with which the whole screen becomes whitish.

When inputting the signal for the horizontal lines

Adjust using the + and - buttons until the lines become horizontal.

When inputting the signal with which the whole screen becomes whitish

Adjust using the + and - buttons until the white color becomes as uniform as possible.

To reset the setting to standard (00)

Press the ↑/+ and ↓/- buttons together.

Troubleshooting

This section may help you isolate the problem. Should the problem persist, unplug the unit and contact your Sony dealer or local authorized Sony service facility.

| Symptom | | Check these items |
|--|--|--|
| No picture | If neither the POWER indicator nor POWER SAVING indicator is lit | Check that the power cord is properly connected. Check that the power switch is set to ON. |
| | If the POWER SAVING indicator is lit | <ul style="list-style-type: none"> • Check the power switch on the video equipment connected to the monitor is set to ON. • Check the input signal and select the corresponding input channel pressing the LINE A/B or RGB/COMPONENT A/B select buttons. |
| If the POWER SAVING Indicator blinks without the picture displayed | | When the monitor detects the abnormal state, the monitor automatically stops operation and makes the POWER SAVING indicator blinks. Turn the monitor power off once, and turn the power on again after about five minutes. If the POWER SAVING indicator still blinks, unplug the monitor and contact your Sony dealer or local authorized Sony service facility. |
| If on-screen menus do not appear | | REMOTE (FRONT LOCK) mode is selected on the REMOTE (RS-232C) menu. For details, see "To control the monitor via the RS-232C interface" on page 22. |
| If the color display is non-uniform | | Trip the power switch once to activate auto-degauss cycle. This function is to demagnetize the CRT to obtain a neutral field for uniform color reproduction. If a seconds degauss cycle is needed, allow a minimum interval of 20 minutes for the best result. Or adjust the landing beam on the LANDING menu. For details, see "To adjust the landing" on page 27. |
| If white does not look white | | Perform color temperature adjustment. For details, see "To adjust the white balance" on page 26. |
| If the screen is too bright and the contrast in the light area is weakened | | Terminate the unit to which the cable connected to the loop-through output connectors on the monitor is connected. This monitor is 75-ohm terminated. However, when you connect the cable to the loop-through output connectors and do not connect the other side of the cable, loop-through output connector on the monitor is unterminated. |
| If picture bounces or has wavy oscillations | | Select the appropriate sync signal on the RGB SYSTEM menu. For details, see "To select the signal to be monitored through the RGB/COMPONENT connectors" on page 22. |
| If the picture is greenish or pinkish (magenta) | | When the picture is greenish: RGB is designated on the RGB SYSTEM menu. Select COMP on the RGB SYSTEM menu. When the picture is pinkish (magenta): COMP is designated on the RGB SYSTEM menu. Select RGB on the RGB SYSTEM menu. For details, see "To select the signal to be monitored through the RGB/COMPONENT connectors" on page 22. |
| If the picture is tilted | | Adjust the landing on the LANDING menu. For details, see "To adjust the landing" on page 27. |

Specifications

Video signal

| | |
|----------------------|--|
| Color system | PAL (Display system: 625/100/2:1) |
| System ^{a)} | 625/50/2:1 (Display system: 625/100/2:1) |
| | 625/50/1:1 (Display system: 625/50/1:1) |

a) System specifications of each system are explained on page 31.

| | |
|------------|--------------------------|
| Resolution | 625/50/2:1: 550 TV lines |
| | 625/50/1:1: 800 TV lines |

| | |
|---------------------|--------------|
| Aperture correction | 0 dB to 6 dB |
|---------------------|--------------|

| | |
|--------------------|-----------------------------|
| Frequency response | 625/50/2:1: 6.8 MHz (−3dB) |
| | 625/50/1:1: 20.0 MHz (−3dB) |

| | |
|-----------------|-------------------|
| Synchronization | AFC time constant |
| | 1.0 msec. (FAST), |
| | 2.0msec. (SLOW) |

Picture performance

| | |
|-------------------------|--|
| Normal scan | 7 % over scan of CRT effective screen area |
| Under scan | 5 % underscan of CRT effective screen area |
| H. linearity | Less than 5.0 % (typical) |
| V. linearity | Less than 5.0 % (typical) |
| Convergence | Central area: 0.5 mm (typical) |
| | Peripheral area: 0.6 mm (typical) |
| Raster size stability | H: 0.5 %, V: 0.5 % |
| High voltage regulation | 0.5 % |

CRT

| | |
|-------------------|-----------------------------------|
| CRT type | HR Trinitron |
| | Aperture grill (AG) pitch: 0.3 mm |
| | EBU standard phosphor |
| Color temperature | 6500 K/5600 K/USER 1/USER 2 |

Inputs

| | |
|----------|--|
| LINE A | |
| VIDEO IN | BNC connector (×1), 1Vp-p +3 dB, −6 dB, sync negative |
| AUDIO IN | Phono jack (×1), −5 dBu ^{b)} , more than 47 kilo-ohms |

LINE B

| | |
|--------|--|
| Y/C IN | 4-pin mini-DIN(×1) |
| | See the pin assignment on the next page. |

| | |
|----------|--|
| AUDIO IN | Phono jack (×1), −5 dBu ^{b)} , more than 47 kilo-ohms |
|----------|--|

RGB/COMPONENT A

| | |
|----------------|-----------------------------------|
| R/G/B channels | BNC connector (×3) |
| | 0.7 Vp-p +3 dB, −6 dB |
| | Sync on green: 0.3 Vp-p, negative |

| | |
|-------------------|---|
| R-Y, B-Y channels | BNC connector (×2) |
| | 0.7 Vp-p +3 dB, −6 dB |
| Y channel | BNC connector (×1) |
| | 1.0 Vp-p +3 dB, −6 dB (standard color bar signal of 75 % chrominance) |

| | |
|----------|--|
| AUDIO IN | Phono jack (×1), −5 dBu ^{b)} , more than 47 kilo-ohms |
|----------|--|

| | |
|----------|--|
| EXT SYNC | BNC connector (×1) |
| | Composite sync 4 Vp-p ±6 dB, sync negative |

RGB/COMPONENT B

| | |
|----------------|-----------------------------------|
| R/G/B channels | BNC connector (×3) |
| | 0.7 Vp-p +3 dB, −6 dB |
| | Sync on green: 0.3 Vp-p, negative |

| | |
|-------------------|---|
| R-Y, B-Y channels | BNC connector (×2) |
| | 0.7 Vp-p +3 dB, −6 dB |
| Y channel | BNC connector (×1) |
| | 1.0 Vp-p +3 dB, −6 dB (standard color bar signal of 75 % chrominance) |

| | |
|----------|--|
| AUDIO IN | Phono jack (×1), −5 dBu ^{b)} , more than 47 kilo-ohms |
|----------|--|

| | |
|----------|--|
| EXT SYNC | BNC connector (×1) |
| | Composite sync 4 Vp-p ±6 dB, sync negative |

| | |
|---------|------------------------------------|
| RS-232C | 9-pin D-sub (×1) |
| | See the pin assignment on page 31. |

b) 0 dBu = 0.775 Vr.m.s.

Specifications

Outputs

- LINE A
- VIDEO OUT
- BNC connector (×1), loop-through,
automatic 75 ohms termination
- AUDIO OUT
- Phono jack (×1), loop-through
- LINE B
- Y/C OUT
- 4-pin mini-DIN (×1), loop-through,
automatic 75 ohms termination
- AUDIO OUT
- Phono jack (×1), loop-through
- RGB/COMPONENT A
- R/R-Y/ G/Y B/B-Y OUT
- BNC connector (×3), loop-through,
automatic 75 ohms termination
- AUDIO OUT
- Phono jack (×1), loop-through
- EXT SYNC
- BNC connector (×1), loop-through,
automatic 75 ohms termination
- Speaker output
- Output level: 0.8 W


General

- Classification of equipment
- Evaluated to EN60601-1, EN60601-1-2
- Type of protection against electric shock:
- Class I equipment
- Degree of protection against harmful ingress of water:
- Ordinary equipment
- Degree of safety of application in the presence of a flammable anaesthetic mixture:
- Not protected equipment
- Mode of operation:
- Continuous operation
- Information concerning type and frequency of technical maintenance:
- Not need maintenance equipment
- Main power switch:
- Functional switch
- Power consumption
- AC 150 W (0.8 to 0.6 A)
- Power requirements
- 220 to 240 V AC, 50/60Hz
- Operating temperature
- 0 to +35°C (32 to 95°F)
- Transport and Storage conditions
- 10 to +40°C (14 to 104°F) Temperature
- 700 to 1,060 hPa Pressure
- 0 to 90% Humidity

- Dimensions
- Approx. 450 ×457.5 × 503 mm
(w/h/d)
(17³/₄ × 18¹/₈ × 19⁷/₈ inches)
not incl. projecting parts and controls
- Mass
- Approx. 34 kg (74lb 15 oz)
- Accessory supplied
- AC power cord (1)
AC plug holder (1)
Side covers (2)
Control panel cover (1)
Panel hinges (2)
Instructions for Use (1)
Interface Manual for Programmers (1)
Quick Reference Card (1)
Double-sided adhesive tapes (4)
Sales Companies Guide (1)


Pin assignment

Y/C IN connector (4-pin mini-DIN)
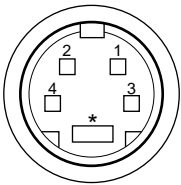
| Pin No. | Signal | Description |
|---------|-------------------------|---|
| 1 | Y-input | 1 Vp-p, sync negative, 75 ohms |
| 2 | CHROMA subcarrier-input | 300m Vp-p, burst Delay time between Y and C: within 0 ± 100 nsec., 75 ohms |
| 3 | GND for Y-input | GND |
| 4 | GND for CHROMA-input | GND |

RS-232C connector (D-sub 9-pin)

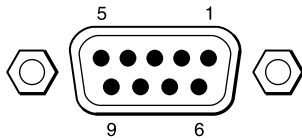
| Pin No. | Signal |
|---------|--------------------------------|
| 1 | — |
| 2 | RX (Remote computer → Monitor) |
| 3 | TX (Monitor → Remote computer) |
| 4 | — |
| 5 | GND |
| 6 | — |
| 7 | — |
| 8 | — |
| 9 | TALLY ON/OFF ^{c)} |

c) ON when pin 5 and pin 9 are shorted.


System specification

Signal timing chart for 625/50/1:1

H (Horizontal)

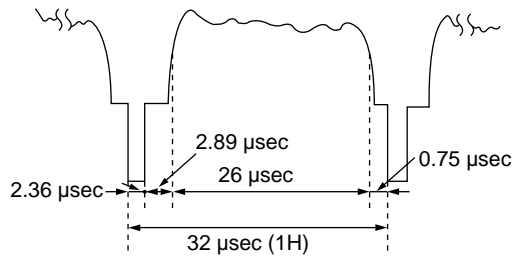
V (Vertical)

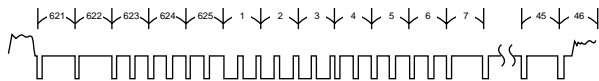
Signal timing chart for 625/100/2:1


H (Horizontal)

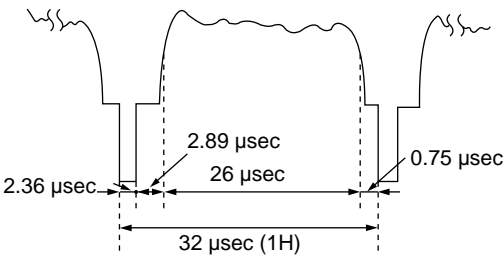V (odd field)

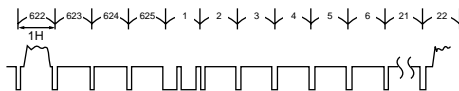V (even filed)

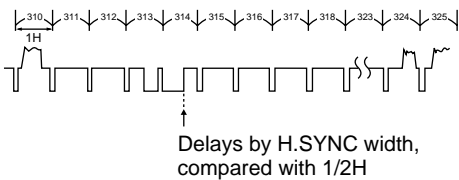Design and specifications are subject to change without notice.

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